



Technical Data Sheet

Thermo Scientific™ Ultra-Low Temperature Freezer

Upright Model Release - 81

Thermo Fisher Scientific, Asheville, North Carolina

Model Number

TDE60086FV

Application, Rating and Electrical Data

Application	Storage of General (non-flammable) Laboratory Materials
Storage Volume	28.8 cu. ft. (816 liters), 600 Standard 2" Boxes
Temperature Rating	-50°C to -86°C
Electrical Power	230V, 50 Hz, 1 Phase
Instrument Rated Current	6.1 AMP
Building Supply Rating	10.0A (min) dedicated grounded circuit
Power Plug / Power Cord	Country Dependant plug / IEC Cords, 10 ft (3.05 m)
Agency Listings	UL, cUL, CE
Application Environment	Non-Corrosive, Non-Flammable, Non-Explosive
	Indoor Use Only, Ventilated 15° C - 32° C (59° F - 90° F)

Dimensions and Construction

Interior Dimensions (H x D x W)	51.2 x 28.3 x 34.4 in. (1300 x 719 x 874 mm)
Exterior Dimensions (H x D x W)	78 x 38.5 x 39.6 in. (1981 x 978 x 1006 mm)
Shipping Dimensions	83.12 x 42.75 x 47.48 in. (2111 x 1086 x 1206 mm)
Shipping Weight / Net Weight	820 lbs. (372 kg) / 709 lbs. (322 kg)
Insulation	Vacuum Insulation Panels with High-Density Water-Blown Polyurethane Foam
Door Seal	Silicone-Based Gasket Seal with Electrical Door Perimeter Heater
Shelves	3 Stainless Steel Adjustable Shelves in 1" (25mm) Increments
Shelf Capacity	Max. Capacity per Shelf: 245 lbs. (110 kg)
Interior / Exterior Material	Painted Steel (Stainless Steel Option) / Painted Steel
All-Direction Casters	Standard with Locks
Other Options	LN2 or CO2 Back Up System, Chart Recorder

Electrical System Configuration



Controller Level	Top
Power Switch	(Rear) Main Circuit Breaker
Controller Type	HIC Interface: 7-segment Display with Capacitive Touch Screen Buttons
Setpoint Security	Yes
Compressor Safe Guard	High Temp Cutout Switch, Current, Logic protection
Control Sensor	Single RTD (1000 ohm Platinum RTD)
Connectivity / Remote Outputs	RS485/4-20mA output/Dry Contacts
Thermo Fisher Cloud	InstrumentConnect™ Remote Monitoring (compatible)
Adjustable Warm/Cold Alarms	Fully Adjustable

Refrigeration Configuration

Refrigeration System	Two Stage Cascade System
Compressor/Number	Industrial Rated, Hermetically Sealed / 2
Compressor Capacity*	559 W
Condenser Type	Enhanced Micro-Channel and Forced-Air Cooled
Expansion Device	Capillary Tube
Evaporator Type	Enhanced Cold Wall Design
Defrost Method	Manual Defrost
Refrigerant (1st/2nd Stage)	R290 / R170+R290 Mix
Environmental Effects	GWP: 3 (R290) / 6 (R170)
	ODP: 0 (R290); 0 (R170)
Flammable	Yes

Typical Performance Characteristics in 20 °C Ambient

Energy Consumption (kW-hr/day)	14.6
Heat Rejection Rate (Btu/hr)	2080
Peak Variation from Setpoint (°C)	+2.8 / -2.4
Sound (dBA)	50.5
1-min Door Opening Recovery to -75°C (min)	24
Average Uniformity at -80°C (°C)	3.6
Average Stability at -80°C (°C)	1.5
Pull Down Time (to -80°C) (hrs)	6.2
Warm Up Time (-80°C to -50°C) (min)	300

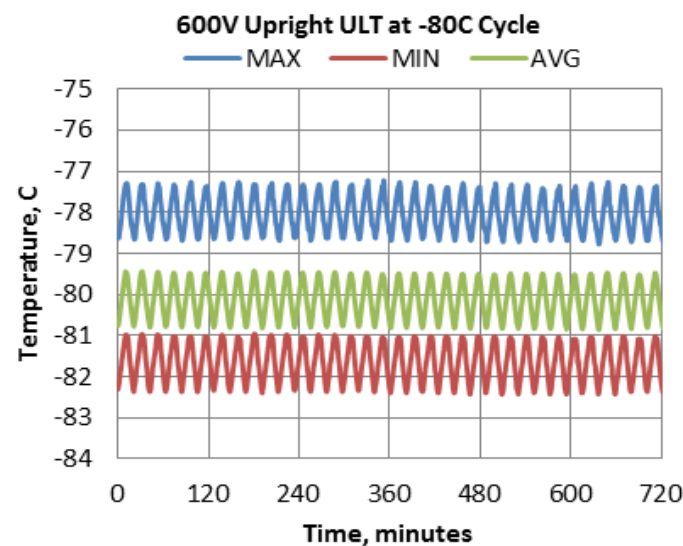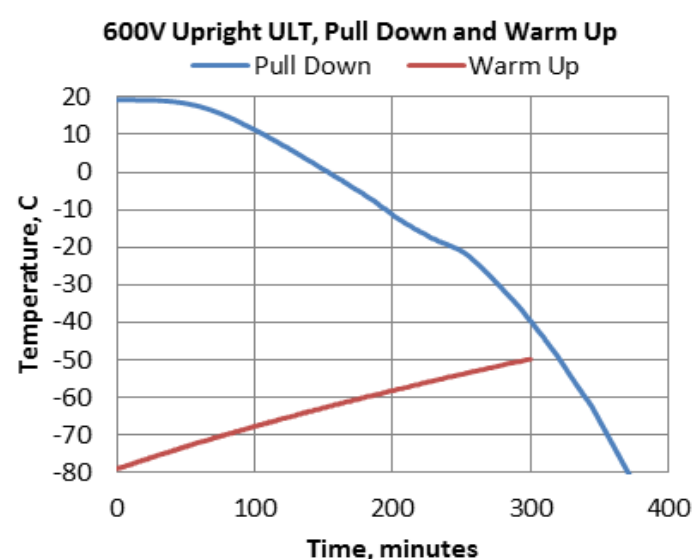

- 1) Performance is nominal and individual units may vary.
- 2) Freezer performance will differ due to product amount, product size and operating conditions.
- 3) Continuous product enhancements may, without notice, result in amendments or omissions to this specification. Thermo Scientific cannot accept responsibility for damage, injury, loss or expenses resulting from misapplication of the information herein.

* Manufacturer measured compressor capacity taken at standard -35°C/45°C (Evap/Cond) condition.

© 2019 Thermo Fisher Scientific Inc. All rights reserved. All trademarks are the property of Thermo Fisher Scientific Inc. and its subsidiaries. Specifications, terms and pricing are subject to change.


Not all products are available in all countries. Please consult your local sales representative for details.



Typical Cabinet Temperature Map
ULT 600, 3 Inner-Shelves + Base, Single Outer Door


Temperatures are averages during > 12
cycles after reaching a setpoint of -80C

Top View of Shelves

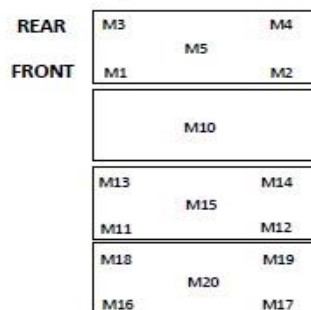

Cabinet Average: -80.1 C
Probe Average: -80.0 C
Peak Variation: +2.8 C / -2.4 C

Front View



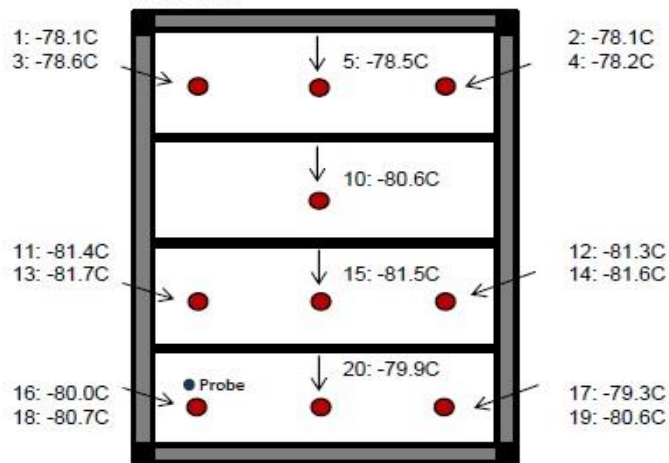


19554-TS-T1


	M1	M2	M3	M4	M5	M10	M11	M12	M13
Avg	-78.1	-78.1	-78.6	-78.2	-78.5	-80.6	-81.4	-81.3	-81.7
Max	-77.4	-77.3	-77.6	-77.2	-77.6	-79.8	-80.7	-80.7	-81
Min	-78.8	-78.9	-79.5	-78.9	-79.2	-81.4	-82.1	-82	-82.4

	M14	M15	M16	M17	M18	M19	M20
Avg	-81.6	-81.5	-80	-79.3	-80.7	-80.6	-79.9
Max	-80.9	-80.9	-79.5	-78.7	-79.9	-79.9	-79.4
Min	-82.3	-82.2	-80.8	-80.1	-81.5	-81.5	-80.5



Thermo Fisher Scientific Proprietary and Confidential

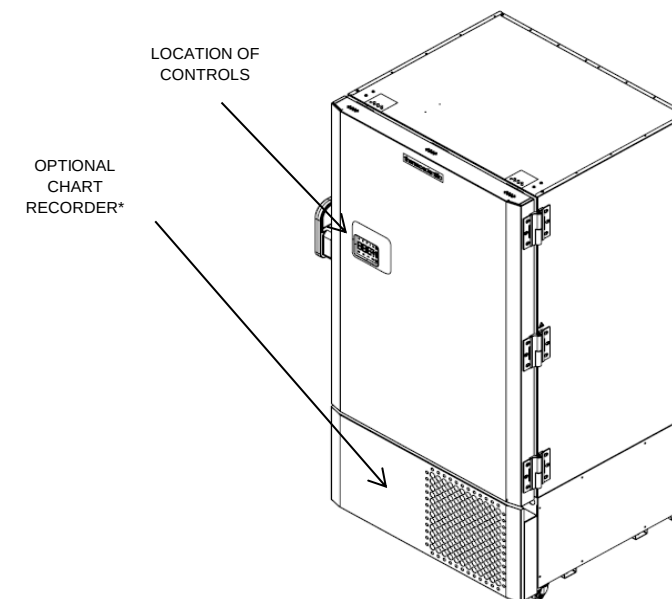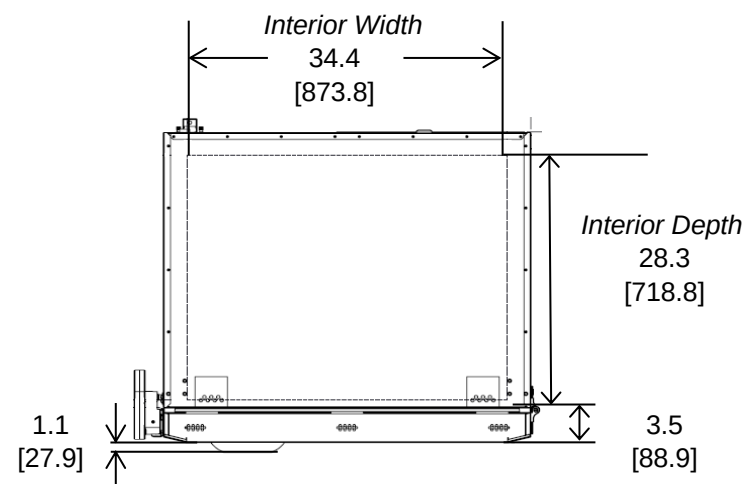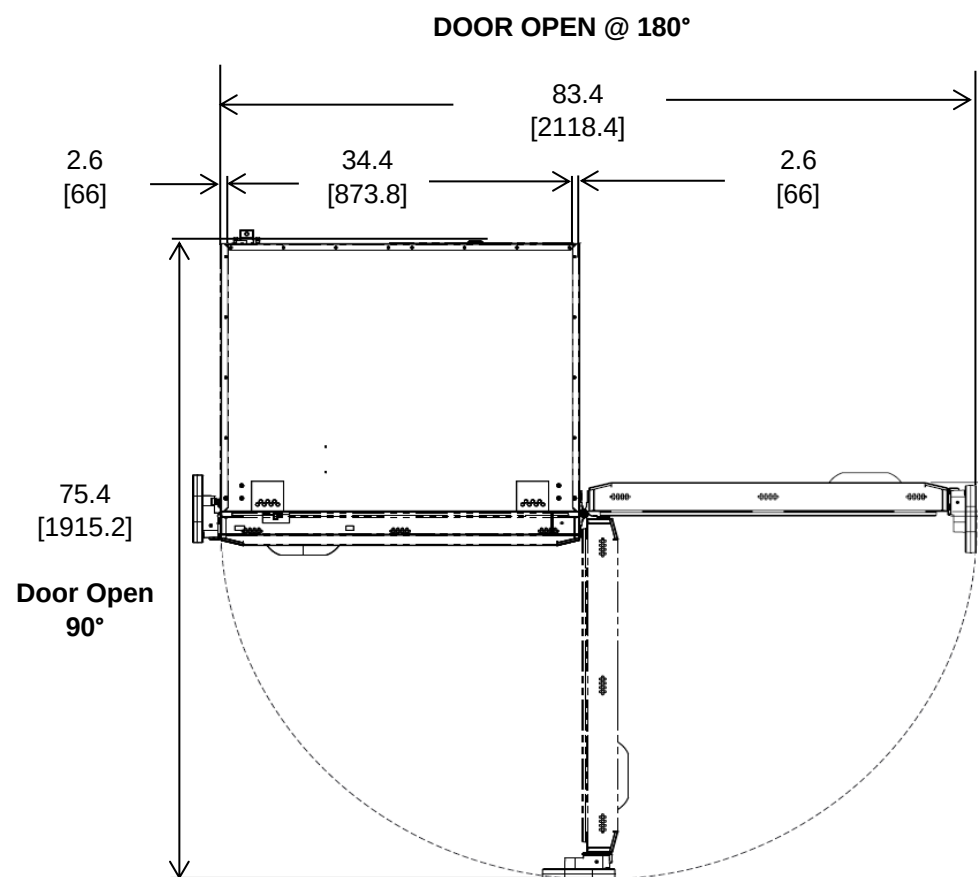

*Chart recorders are stand alone on 300/400 models
*Backup systems are stand alone on 300 models.

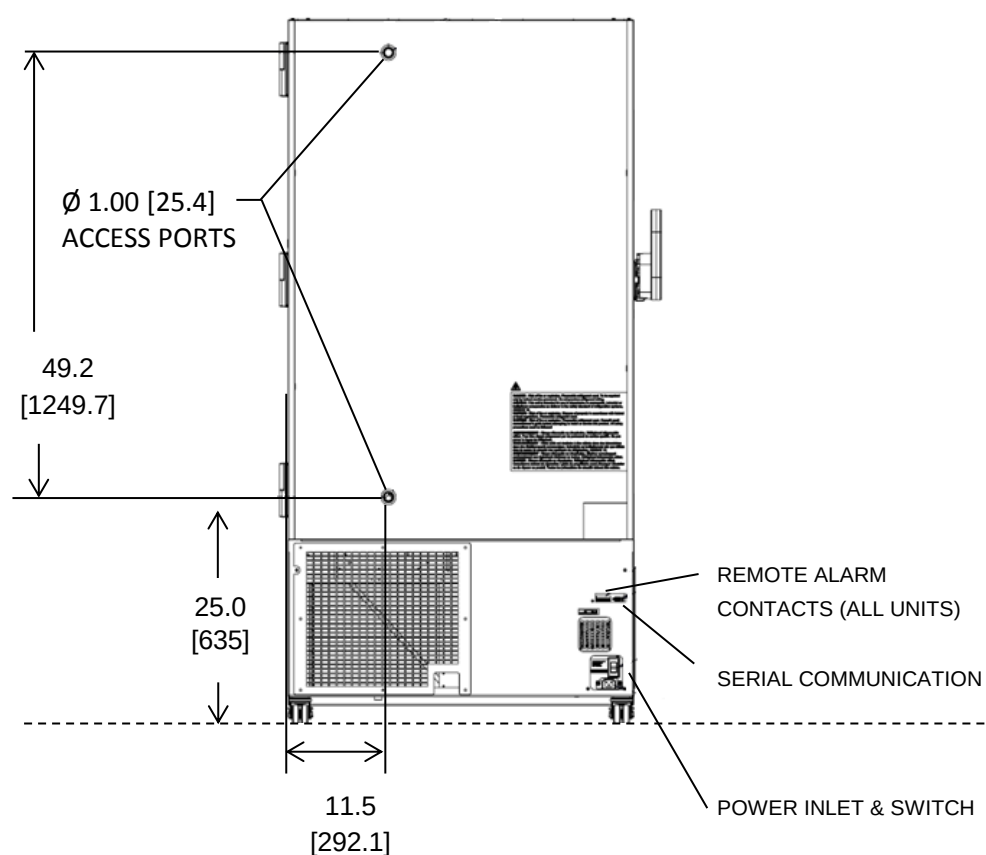

BACK

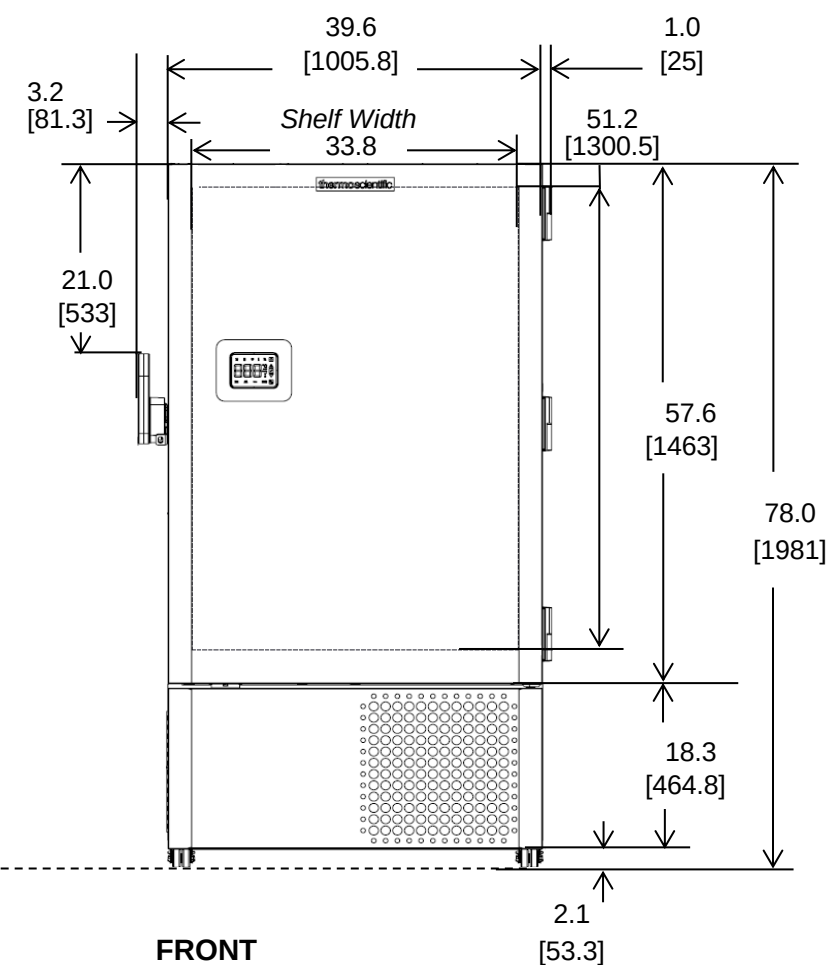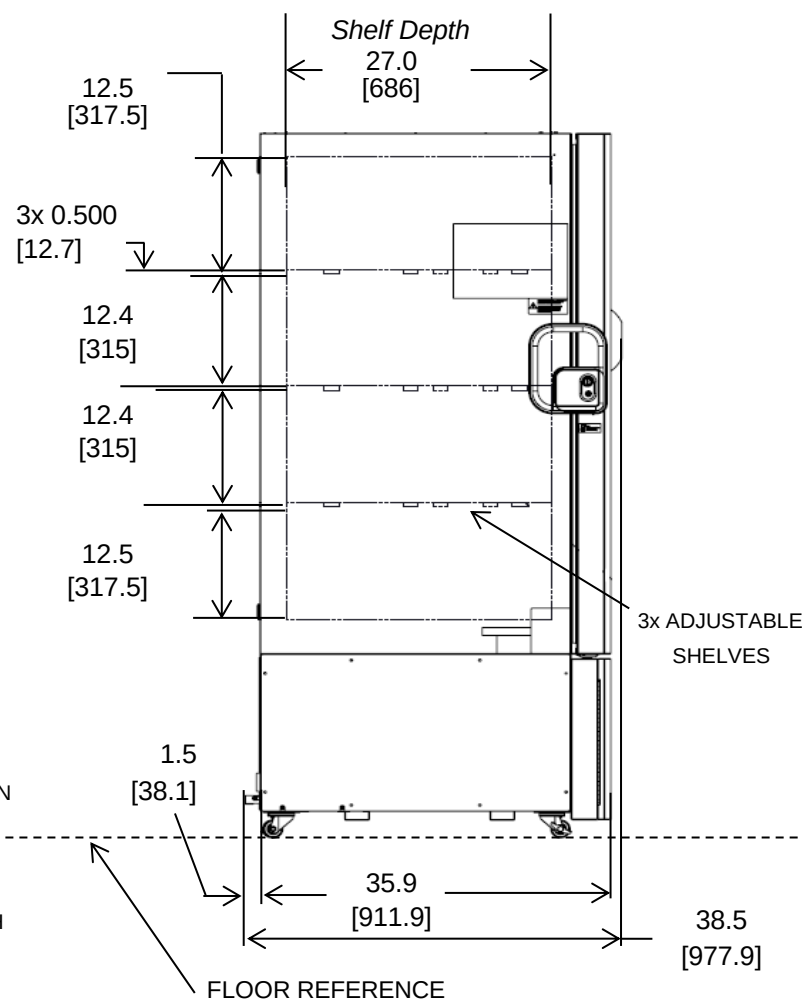

NOTE: DUAL DIMENSION IS INCH OVER METRIC

DO NOT USE FOR ENGINEERING PURPOSES. SUBJECT TO CHANGE WITHOUT NOTICE.

thermo
scientific

28.8 cuft Upright Freezer
4 Inner Doors
Single Outer Door
Top Mount Controls